


# MDO-C Alveolar Distractor

Simple and effective

Developed in cooperation with Albino Triaca, MD, DMD and Roger Minoretti, MD, DMD,  
Center for Maxillofacial Surgery, Pyramide Clinic, Zürich, Switzerland

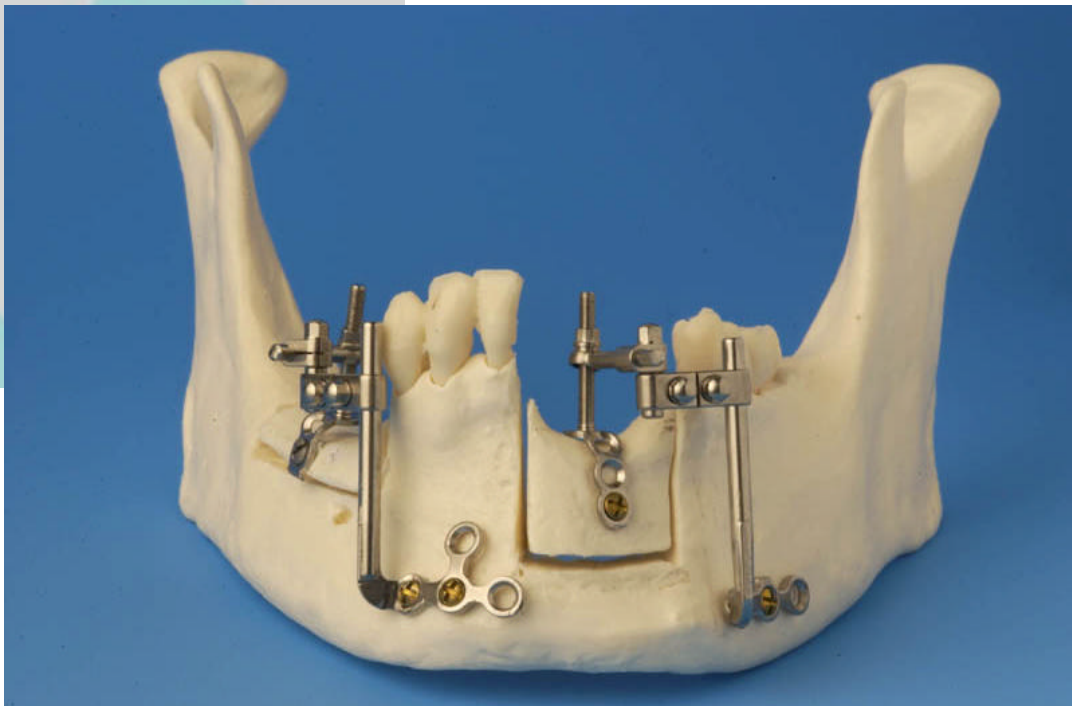

- Remote fixation, helps to avoid nerve lesions even in very atrophic cases
- Extra-mucosal placement of mechanics, simple surgery, less scar tissue.
- Correction of vector during distraction. Prevention or correction of lingual deflection.
- Helps to create bone where it is needed

Patent pending

## Orthognathics MDO-C – Alveolar Distractor

The Orthognathics MDO-C Alveolar Distractor allows a versatile application for pre-prosthetic segmental distraction in the mandible or the maxilla.



The MDO-C Alveolar Distractor is a simple yet very effective tool for the distraction of alveolar defects. Due to its modular approach, it allows optimal adaptation to the specific situation as well as easy installation and confi-

guration. The possibility to modify the distraction vector during activation allows to avoid unwanted tipping towards lingual.

The fixation remote from the distraction zone helps to avoid any problems with anatomical structure such as nerves.

### Technical Data:

Distraction range	15 mm
Material	surgical steel
Bone plates	for 1.5/2 mm screws

For any questions on this product, please send us an e-mail or write us by fax or mail.

- Send more information
- Send price quotation
- Please arrange a meeting or demonstration/surgery
- Please keep me informed about workshops and presentations on MDO-P and other Orthognathics products

Orthognathics Ltd.  
 Winterthurerstr. 92  
 CH-8006 Zurich/Switzerland  
 Tel. +41 (0)43 2689200  
 Fax +41 (0)43 2689202  
 info@orthognathics.com

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Phone: \_\_\_\_\_

Fax: \_\_\_\_\_

e-mail: \_\_\_\_\_

\_\_\_\_\_  
 Date Signature